

USER'S GUIDE



Your Oxygen Facial System includes:

- Main unit
- Removable and adjustable mirror
- 3 fl oz oxygenating foaming cleansing milk
- 3 fl oz oxygenating moisturising cream
- Separate recharger adapter
- Soft carry case



THE NEW
OXYGEN FACIAL SYSTEM
FROM VERSEO

THE NEW OXYGEN FACIAL SYSTEM FROM VERSEO

The easy way to a healthier complexion!

- Dries out blemishes • helps reduce fine lines • soothes puffy eyes •
- decreases open pores • enhances the results of oxygenating skin care products •
- rechargeable so can be used anywhere •

The Verseo Oxygen Facial System will enhance the appearance of your complexion. Oxygen technology has been used worldwide for over ten years.

The skin is continually under attack from the ever deteriorating environment, chemical pollutants, smoke, dust and everyday pollution. This can lead to premature ageing, clogged-up pores and unsightly blemishes.

The Oxygen Pen uses a direct high frequency current (DHF) bulb which is gently applied to the skin. The DHF bulb lights up near or in contact with the skin. This is due to the DHF bulb ionising the oxygen in the air, forming ozone gas enriched with oxygen (O₃). The ozone gas gives off a slight harmless smell.

Enriched oxygen (O₃) dramatically improves the skin condition. oxygen is nature's cleanser creating healthy dermal nutrition, harmful toxins and bacteria find it hard to survive in an oxygen rich environment.

Within a short time your skin will look smoother and your complexion firmer.



WHAT IS THE SECRET OF THE OXYGEN FACIAL SYSTEM

....in a word Oxygen.

The skin needs an ample supply of oxygen to repair and maintain its self, after the age of 20 the skin begins to lose its ability to absorb oxygen.

Studies show that the skin on the face of a 25 year old is as much as 30 - 40% thinner than that of an adolescent or child.

What influences ageing skin

For a start environmental and behavioural practices.

Over exposure to UV radiation, pollution, smoking, flying and direct contact with irritating substances such as cleaning products, and of course stress. Reduced circulation and decreasing capillary action being the most important factors.

What are the results of these ageing influences

Progressive decrease of micro circulation resulting in a lower oxygen supply making the skin cells less active. The results are, lower fat secretion and less sweat so that skin becomes dryer and cell renewal slows down. Both collagen and elastin production decreases.

HOW DOES THE VERSEO OXYGENATING MOISTURISING CREAM WORK

The older we get the partial oxygen pressure in the capillaries decreases. (Particularly on the face because of exposure to the elements and the less extensive capillary network).

Therefore early ageing sets in and the skin becomes thinner with the appearance of lines and wrinkles, age spots and poor complexion. The face sags, as there is less collagen and elastin.

The problem is how to deliver the oxygen and nutrients (such as vitamin A & E which help to nourish the skin and underlying collagen layer) to the skin to supplement the loss.


The answer comes from a sophisticated cream containing hydrogen peroxide. This type of formulation was originally used to help the recovery of major skin burns and post operative plastic surgery. The result helps make the skin more elastic and softer which will enhance your appearance.

The cream contains stabilised oxygen in the form of hydrogen peroxide as well as neutral vitamin A. When applied to the skin the oxygen atom (O) becomes an oxygen molecule (O₂). The molecules aggregate and the liquid oxygen changes to gas. This pressure release has the effect of micro massage on the outer layers of the skin, unblocking capillaries. Likewise the vitamin A becomes acidified helping to stimulate collagen.

The skin firms and strengthens, the colour becomes more even, and the pores reduce in size with fine lines and wrinkles smoothed out.

The combined use of the Verseo Oxygen Facial System which creates ozone (enriched oxygen O₃) with the direct massage effect enhances the performance of all oxygenating products.

PRODUCT GUIDE

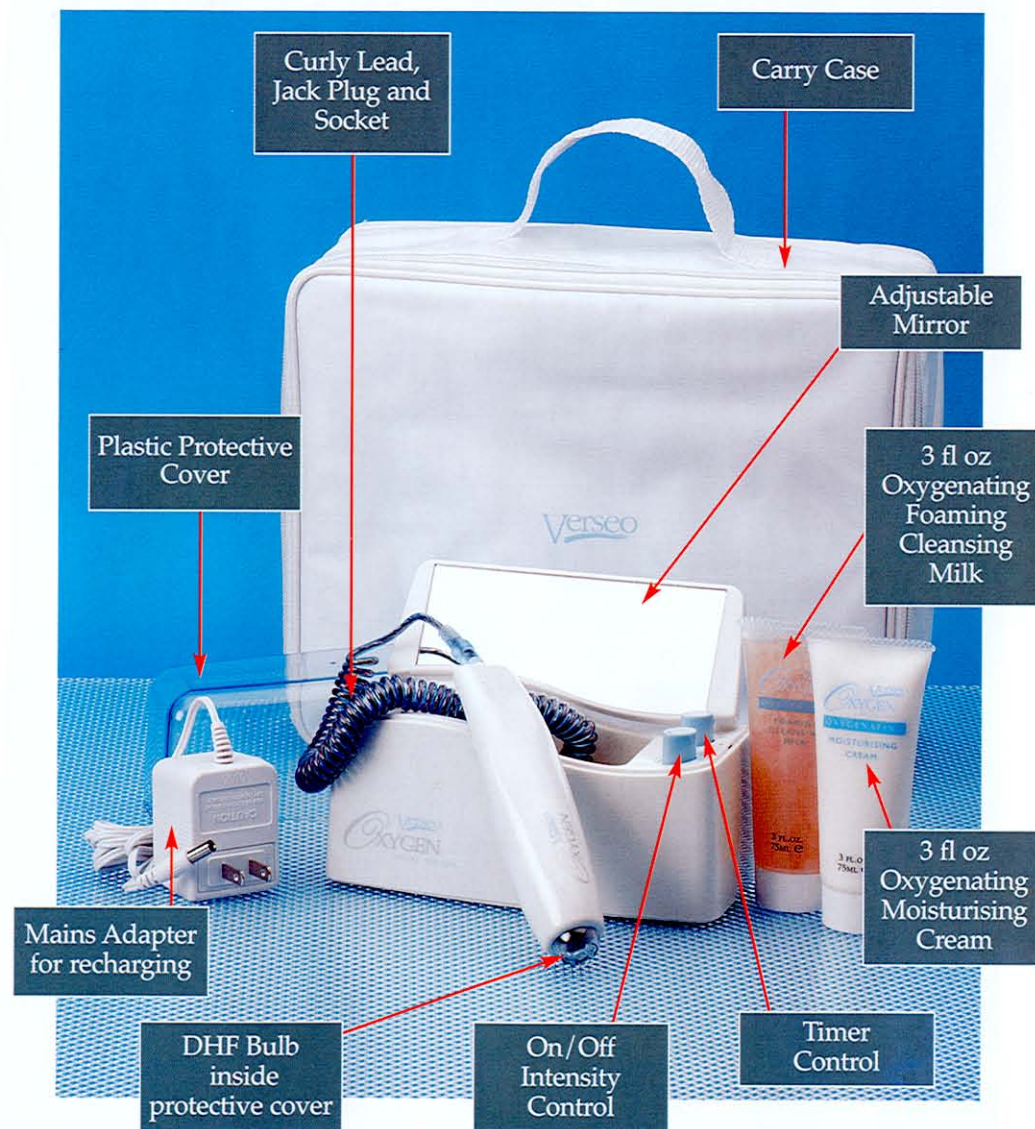

HOW TO USE THE OXYGEN FACIAL SYSTEM

Pre Treatment

1. Remove all traces of dirt or make up by thoroughly cleansing the skin with non alcoholic or PH balanced cleanser or use the Verseo Oxygenating Foaming Cleansing Milk supplied with this product. Follow the instructions on tube.



2. Test the Oxygen wand by gliding the bulb over the back of your hand first. Turn Timer ON, and then turn intensity ON.

Turn the intensity up to its highest setting maintaining comfort. The bulb will crackle and spark. (Absolutely safe). The harmless smell of enriched oxygen (ozone) will be present after contact with the skin. Follow this procedure when applying the wand to your face.



3. When concentrating on blemishes it sometimes helps to lift the Oxygen wand just clear above the skin surface. This allows a more focused stream of energy onto the blemish.

Starting Treatment

1. Plug recharger adapter lead into jack socket at the back of unit. (Use only recharger adapter Model No.ILD35-120200 supplied with product).
2. Plug recharger adapter into mains.
3. When charging product for the first time, charge for minimum 6 hours. Charger will automatically cut-out when battery is fully charged. After initial 6 hour charge, battery will only take 2-3 hours to recharge.
4. A fully charged battery should provide approximately 10 treatments of 3 minutes duration.

5. Turn ON/OFF intensity control to 'ON' to lowest setting 1. Leave on until TIMER control turns unit off automatically so DHF bulb is warmed up before treatment.

NOTE: The TIMER control is marked 1 to 5. Each number is marked approximately (but not accurate) equivalent to a minute i.e. 3 = 3 minutes

6. Now start treatment on selected area using:

LOW INTENSITY 1 Minute Treatment

7. Run Oxygen wand on or as close to skin as possible. Turn up intensity to the highest setting without feeling discomfort (the higher the setting, the better the results). Start with a low setting and turn up gradually. See chart for optimum time setting. Use the Oxygen wand in a circular motion.
8. For blemishes, hold Oxygen wand over blemish for 2-3 seconds, a good energy flow produces sparks which will immediately dry out the blemish.
9. Once treatment is complete, thoroughly moisturise with oxygenating cream provided. Follow instructions on tube.
10. The DHF bulb is enclosed underneath protective cover. In the unlikely event bulb fails call our Helpline on 562-948-1114.
11. If you use in bathroom, make sure recharged unit is disconnected by removing jack plug from back of unit. Take care to keep away from water. NEVER USE IN BATH.



GUIDANCE CHART FOR OPTIMUM TIME AND INTENSITY SETTINGS

The chart below is for guidance only. The higher the setting, the better the result. (Always acclimatise the treatment area first on a low setting). Experiment and see which programme suits you best.

If excessive peeling occurs, reduce intensity setting. If irritation or discomfort is felt followed by redness, reduce the time of treatment.

Skin Type	Maximum Duration	Daily Use	Intensity Level
Fair Sensitive	1 minute	1	Low
Fair Oily	2 minutes	1	Medium
Average Dry	2 minutes	1	Medium
Average Oily	3 minutes	1	High
Mature	3 minutes	1	Medium



IMPORTANT INFORMATION BEFORE USING THE OXYGEN FACIAL SYSTEM

Do's & Don'ts:

- If DHF bulb fails return to Atlantic Horizon. We will replace wand free of charge.
- Do not use outdoors
- Do not immerse in water
- Do not use near inflammable or combustible material
- Do not use if pregnant
- Do not use near eyes, mouth or nose
- Do not use acne treatment prior to using the Oxygen Facial System
- Do not use if you are suffering from diabetes, phlebitis, psoriasis or epilepsy
- Do not use if you have any heart condition or use a heart pacemaker
- Do not exceed recommended time for maximum use (see page opposite), or use more than once in 24 hours
- Do not use alcohol based cleansers prior to treatment
- Do not use perfume or perfumed products prior to use
- Do not use on inflamed skin, sunburn, open abrasions or on moles
- Do not use near varicose veins or broken capillaries

If in doubt consult your doctor before use

Your Questions Answered

How frequently can I use the Oxygen Facial System?

Do not use on the same treatment area more than once a day, unless you have oily skin. If you use twice a day on a normal skin the Oxygen Facial System will dry out the skin.

Is the Oxygen Facial System painful to use!

No. During treatment there is a slight tingling sensation. This is normal and completely safe.

Are there any side effects?

The skin may go slightly pink due to the increase of blood flow to the surface of the skin. This disappears shortly after treatment.

When can I apply make up?

Wait for at least one hour after treatment.

Why does my skin peel?

This means the intensity is too high and the time of treatment is too long. Exfoliation is created by the treatment but skin peeling is the result of excessive treatment.

What product do I use to thoroughly cleanse my skin before treatment?

Use the Verseo Oxygenating Foaming Cleansing Milk provided or any PH balanced cleansing product.

What products do I use after treatment?

The Verseo Oxygenating Moisturising Cream, containing Vitamin E and A, Allantoin, Sodium Hyaluronate and Hydrogen Peroxide or any quality Oxygenating Product.

The Oxygen Facial System doesn't work?

Send to customer services for repair/renewal (see address opposite).

IMPORTANT INFORMATION

General Information

Do not immerse the Oxygen Facial System in water.

Technical Data

Model No: 1602US

Input Voltage 120V - 60 Hz

Output Current: 12V 200Ma

Complies to: EN50081 - 1: 1992, EN55022: 1994, EN6100-3-2: 1995, EN500H2-1 1997.

PART 15 OF THE FCC RULES.

Conforms to: EMC Regulations (Electro-Magnetic Compatibility Directive 89/336 EEC)

Use only UL approval recharger adapter supplied by Atlantic Horizon Model No. ILD35-120200.

Repairs and Services

If the unit is damaged, return to Atlantic Horizon. Do not tamper with wires or attempt to open control unit as this will invalidate the guarantee.

Guarantee

The Verseo Oxygen Facial System is guaranteed for one year against original factory faults in material and workmanship. If repairs are necessary as a result of accident or misuse is not included in our Guarantee.

Important: This undertaking is in addition to consumer's statutory rights and does not affect those rights in any way.

The Verseo Oxygen Facial System is manufactured in China to the exclusive USA design of Atlantic Horizon International Inc., 110W 34th Street, Suite 1106, New York NY10001.



TREATMENT RECORD

Keep a record of your treatment times, intensity settings and any notes for your reference when commencing new treatment.

Date: _____ Time: _____ Intensity Setting: _____

Notes: _____

Date: _____ Time: _____ Intensity Setting: _____

Notes: _____

Date: _____ Time: _____ Intensity Setting: _____

Notes: _____

Date: _____ Time: _____ Intensity Setting: _____

Notes: _____

Date: _____ Time: _____ Intensity Setting: _____

Notes: _____

Date: _____ Time: _____ Intensity Setting: _____

Notes: _____

Date: _____ Time: _____ Intensity Setting: _____

Notes: _____

Date: _____ Time: _____ Intensity Setting: _____

Notes: _____

Date: _____ Time: _____ Intensity Setting: _____

Notes: _____

Date: _____ Time: _____ Intensity Setting: _____

Notes: _____



ORDER FORM

Please send me the following items listed below:

3 fl oz Oxygenating Foaming Cleansing Milk	Quantity	Price:\$19.99
Item No.1623		
3 fl ozl Oxygenating Moisturiser Cream	Quantity	Price:\$19.99
Item No. 1624		
Shipping and Handling		\$3.95
		Total Amount:

Name: _____

Address: _____

Postcode: _____

To order by phone: Verseo Customer Service at 562-948-1114 (Credit Cards accepted).

Or Send Check or Money order to: Verseo Inc. PO Box 2727,Downey, CA 90242.

"For our latest offers, you can visit us at www.Verseo.com."

ORDER FORM

Please send me the following items listed below:

3 fl oz Oxygenating Foaming Cleansing Milk	Quantity	Price:\$19.99
Item No.1623		
3 fl ozl Oxygenating Moisturiser Cream	Quantity	Price:\$19.99
Item No. 1624		
Shipping and Handling		\$3.95
		Total Amount:

Name: _____

Address: _____

Postcode: _____

To order by phone: Verseo Customer Service at 562-948-1114 (Credit Cards accepted).

Or Send Check or Money order to: Verseo Inc. PO Box 2727,Downey, CA 90242.

"For our latest offers, you can visit us at www.Verseo.com."

ORDER FORM

Please send me the following items listed below:

3 fl oz Oxygenating Foaming Cleansing Milk	Quantity	Price:\$19.99
Item No.1623		
3 fl ozl Oxygenating Moisturiser Cream	Quantity	Price:\$19.99
Item No. 1624		
Shipping and Handling		\$3.95
		Total Amount:

Name: _____

Address: _____

Postcode: _____

To order by phone: Verseo Customer Service at 562-948-1114 (Credit Cards accepted).

Or Send Check or Money order to: Verseo Inc. PO Box 2727,Downey, CA 90242.

"For our latest offers, you can visit us at www.Verseo.com."